
1a


1b

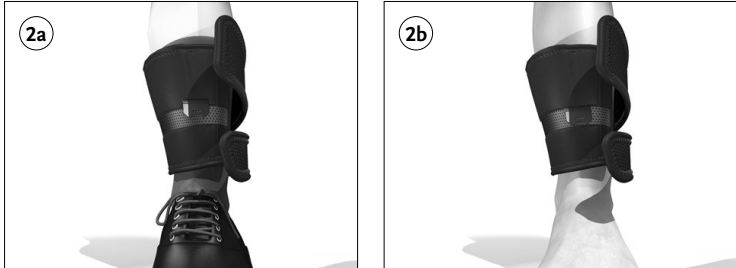2a

2b

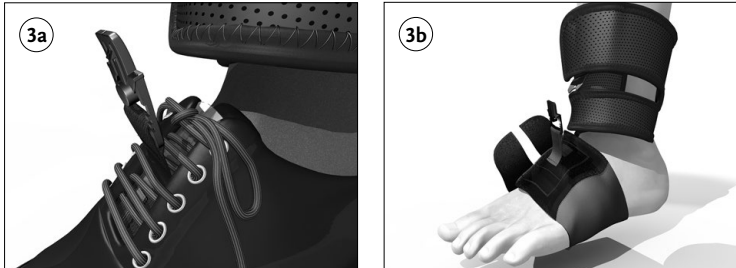3a

3b

4

ENGLISH

MD Medical Device

INTENDED USE

The device is intended to support drop-foot.

Indications for use

Drop foot.

No known contraindications.

Warnings and Cautions:

If you have a medical condition that decreases blood circulation in your limbs, then consult your healthcare professional before you use the device. Such medical conditions include, but are not limited to, diabetes or peripheral vascular disease.

GENERAL SAFETY INSTRUCTIONS

Read these instructions carefully before use. Keep them for future reference. Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities. The patient should stop using the device and contact a healthcare professional:

- If there is a change or loss in device functionality, or if the device shows signs of damage or wear hindering its normal functions.
- If any pain, skin irritation, or unusual reaction occurs with the use of the device. The device is for single patient – multiple use.

FITTING INSTRUCTIONS

Device Application
There are two versions, with the inlay (A) for in-shoe use and the foot wrap (B).

- Measure the circumference of the ankle around the malleoli (Fig. 1a).
- Detach the inlay from the ankle strap by pressing the button in the middle of the quick-acting closure.

DISPOSAL

The device and packaging must be disposed of in accordance with respective local or national environmental regulations.

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